

# Zeencare

## HEATED SHOULDER WRAP USER MANUAL












## FEATURES



1. Heated Shoulder Wrap
2. Timer Setting
3. Heat Setting
4. USB Type C connection port
5. USB Cable
6. Powerbank Pouch
7. Adjustable Strap

# HOW TO START

- 1 Power On and Power Off**  
To turn on , just plug the heated wrap to power bank or other USB power supply. To turn off, just unplug it.
- 2 Switching Temperature**   
After powering on, press the heat button  for 3 seconds to switch to high level temperature , press again for medium level and  press again for low level. 
- 3 Timer Set up**   
After level temperature is chosen, Press the timer button  once for 120 mins. Press the  again once for 60 mins and press the  again once for 30 mins.

## **IMPORTANT:**

This product is not rechargeable and requires an external power source. If it does not turn on, we recommend switching to a different power supply as some may not be compatible.


# PRODUCT SPECS

Net Weight : 180 g

Battery Power: 2A, 9W

Rated Voltage: 5.0V

Power Socket: Type-C

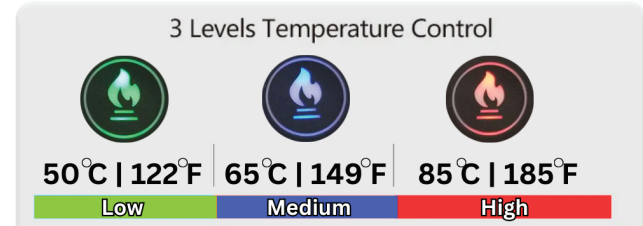

## Precautions for Use

Don't Hand Wash	Don't Machine Wash	Do Not Rip	Do Not Dry Clean	Do not Bleach
